HOMEDICS®

MiBody 360 Personal Scale

Quick Start Guide

360SC-EU

12 year
guarantee

Quick Start Guide

TO DOWNLOAD FULL INSTRUCTION MANUAL AND INSTALL MIBODY 360 TRACKING SOFTWARE

PLEASE VISIT:

www.mibody360.com

TO START USING YOUR HOMEDICS MIBODY360 PERSONAL SCALES STRAIGHT AWAY SIMPLY FOLLOW THE QUICK START INSTRUCTIONS BELOW:

IMPORTANT: DO NOT IMMERSE THE SCALE IN WATER.

The scale can store up to 30 results per person. Results are transferred to your computer using the USB flash drive supplied. Once installed you can analyse and manage your results using the MiBody360 computer software.

Getting started:

Step 1: Batteries

Insert the 4 x AA batteries supplied into the battery compartment, ensuring the + and – terminals are correctly orientated.

Step 2: Carpet feet

For use on carpet please attach the carpet feet provided (UK only).

User set up:

TO SET UP EACH INDIVIDUAL USER YOU WILL FIRST NEED TO DOWNLOAD THE MIBODY 360 SOFTWARE FROM:

www.mibody360.com

Once you have set yourself up as user on the software by entering all of your personal details you must set yourself as the **corresponding user** on the scale.

You will then be able to take full advantage of all the software's features and displays including setting goals.

For full instructions on this please download the complete instruction booklet from www.mibody360.com

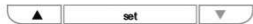
insert batteries press  or  to set day then press 'set'. Repeat for month, then press 'set' and year followed by 'set'

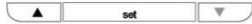
press  or  to fix the hour, then press 'set'. Do the same again for minutes and press 'set'

Hold 'set' for 2 seconds to confirm your settings

press  to enter your user number
 IMPORTANT: select same user number you have allocated yourself on PC software

Step 3: Select Weight Mode.

Note: Not all models have this feature.
 European models standard is set to KG.

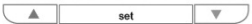
In use:

use  or  to enter user number
 Step on the scale and the display will show your weight
 all other measurements are stored in the memory and transferred to your PC via the USB memory stick


to download to the USB memory stick insert the stick then press 'set'
 data from the scale is uploaded to the USB memory stick when uploading is finished the display will show 'End'
 DO NOT remove the USB stick during uploading

if the memory reaches it's capacity then the display will show 'FULL'
 when a user steps on the scale

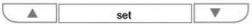when downloaded the display will return to normal

Weight Reading Only – or guest use (note: readings will not be saved to memory for transfer to your computer)

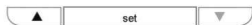
If you just want to know your weight, follow this sequence:

Stand on platform.
 Weight is displayed.

Once measurements have been taken after 20 seconds the scale will revert to the default display of date and time only.

step on scales to read weight in kgs

WARNING MESSAGES:

when the batteries are low the display will show 'bat lo' to highlight the condition

Indicates the batteries are low, simply replace the 4 x AA batteries into the battery compartment, ensuring the + and – terminals are correctly orientated

any other fault will be displayed as an 'error' message on the display

For full instructions on downloading full memory to PC please follow complete instructions booklet from www.mibody360.com

Information for your safety:

This product is for domestic use only.

Professional medical guidance should always be sought before embarking on diet and exercise programmes.

The HoMedics MiBody360 personal scale is calibrated for adult use (age 16+ years). It is not suitable for use during pregnancy.

CAUTION: Do not use this product if you have a pacemaker or other medical device fitted. If in doubt, consult your doctor.

■ This marking indicates that this product should not be disposed with other household wastes throughout the EU. To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources. To return your used device, please use the return and collection systems or contact the retailer where the product was purchased. They can take this product for environmental safe recycling. (WEEE explanation)

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